



Instructions for reassigning a unit to another system


1. After you have logged on select “Change System” from the File Menu

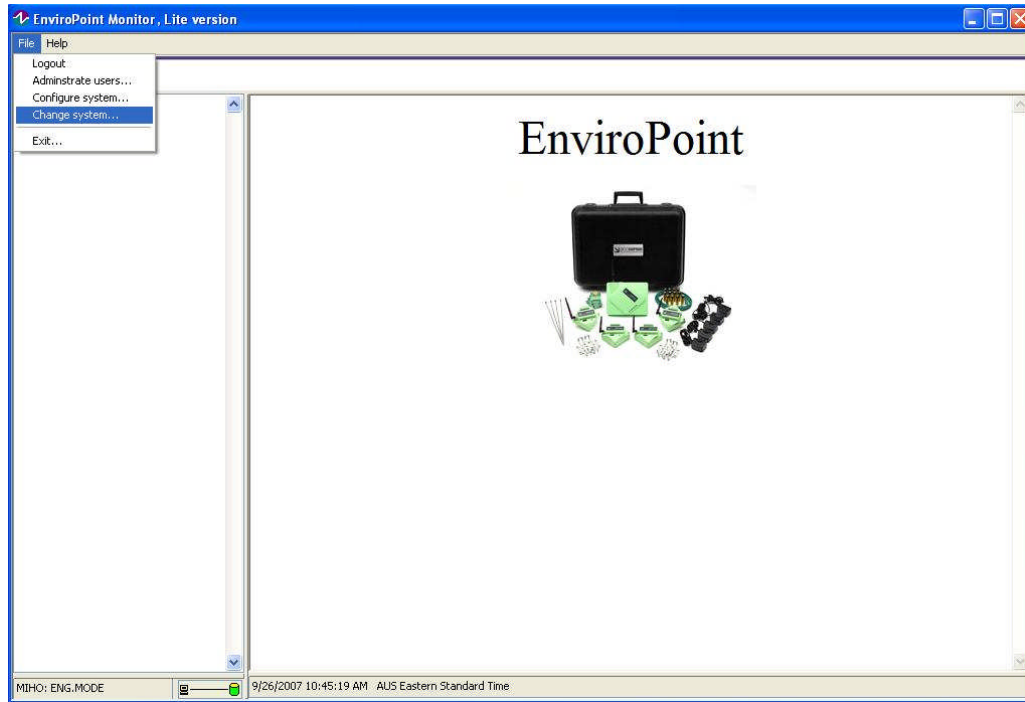


2. In the dialog box that appears select the System “ANMM Trial” and press OK
3. Log in to the new system and then Select “Configure System” From the File menu
4. Expand the menu tree under the system name by clicking the + symbol
5. Select Delete Point From the Shortcut menu (right click) on the point name

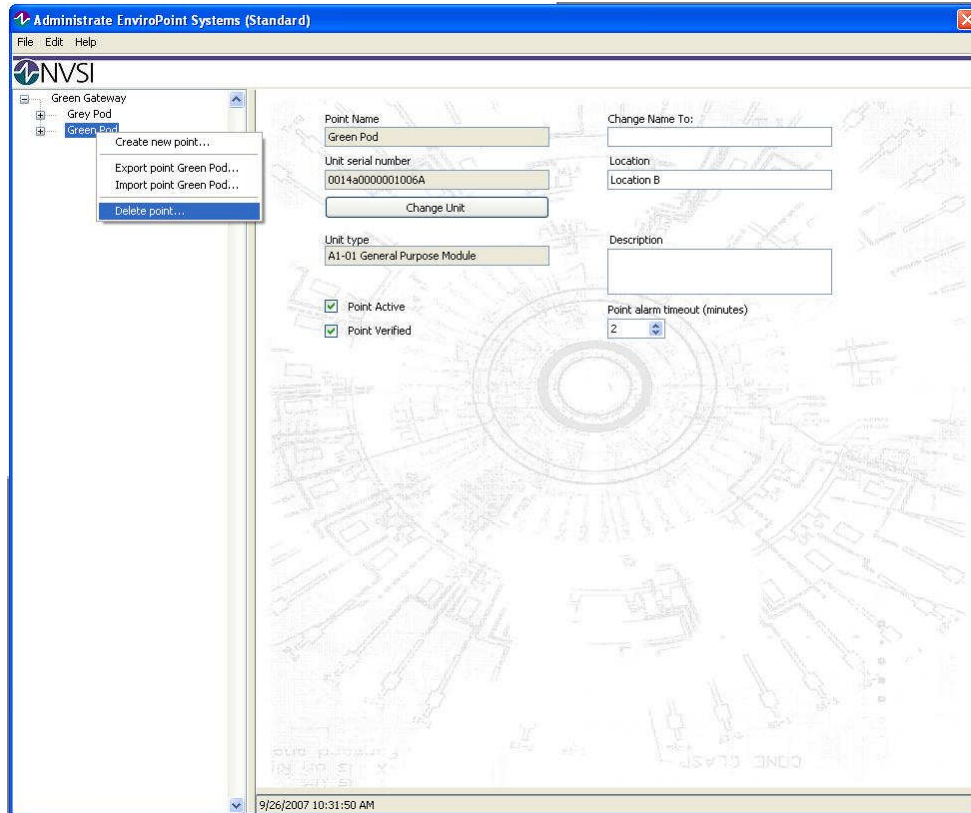


6. Select “Save” From the File menu

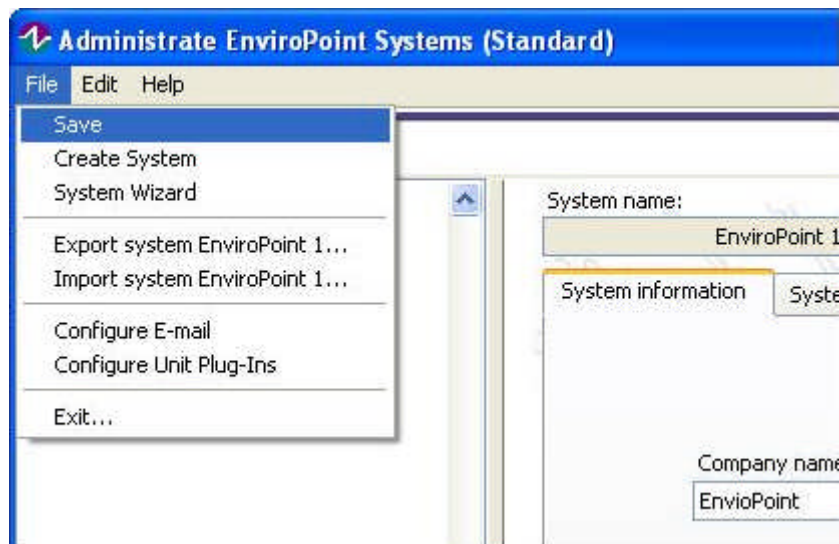


7. Enter change number and description. When save is complete, exit from the configuration screen back to the main program and apply settings when prompted

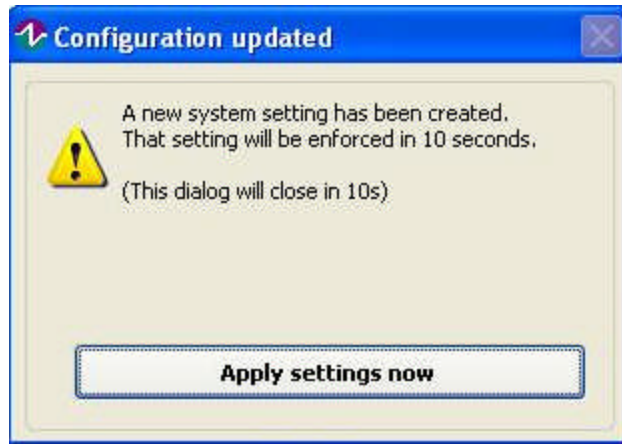


8. Repeat the steps 1-3 to get back into the "ANMM Trial 2" system configuration
9. Select "Create new point..." from the Shortcut menu (right click) on the system name.

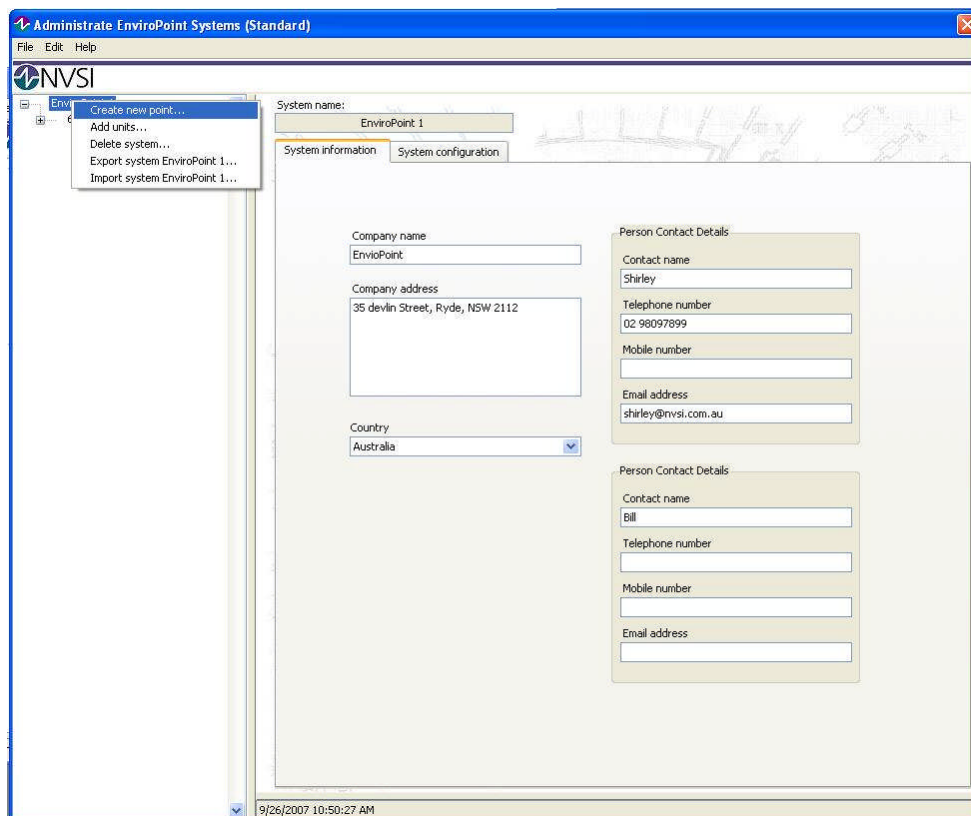

10. Select the unit Serial Number that you wish to add from drop menu, enter an appropriate point name and select OK.

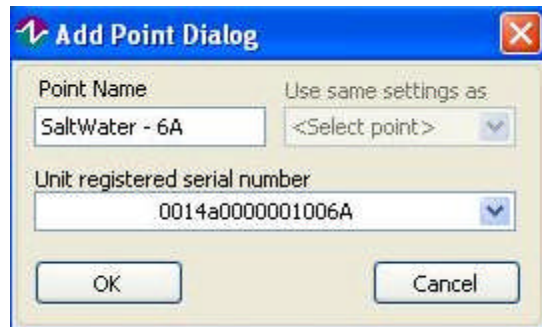

11. Expand the tree to show the new point name, click on the point name and make sure the verified box is ticked on the main panel

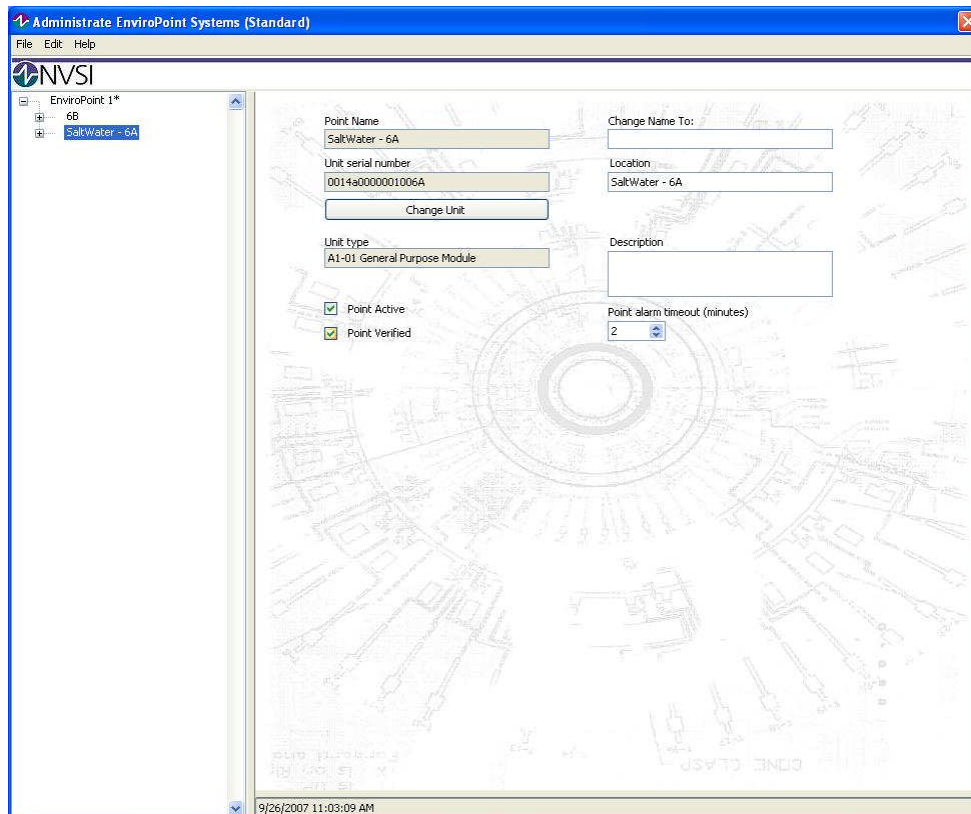

12. Finally save your system as before.
13. After you apply the new settings you should restart the application to make sure all of the data is refreshed properly.

EnviroPoint Save ✕

Save System!

You are about to save the current system settings

This may take up to a few minutes and will cause other monitors to stop displaying live data

Station ID

Change Number:

Change description:

Verified Set Sensor CalibLock to Permanent